
POWER WHITENING KIT

Wireless Accelerator & 6 x Power Whitening Gel

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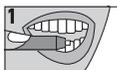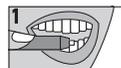
Made in P.R.C.

6 x 5ml e

The Power Whitening Kit ensures maximum bleaching results and up to 70% whiter teeth. The set includes 6 Power Whitening Gel Pens for maximum results in a few days. The wireless Whitening Accelerator with 32 Red and Blue WaveLight LEDs is suitable for use at home or on the road and can be easily charged using the included cable.

Application

1. Clean your teeth with a toothbrush.
2. Press the end of the pen until bleaching gel comes out.
3. Apply the gel evenly to all visible teeth. Make sure that the gel does not touch your gums.
4. Enclose the whitening LED with your lips and start the teeth whitening process with the timer button.
5. Depending on the degree of discolouration, you can let the bleaching gel work for 15–30 minutes.
6. Remove the accelerator and rinse your mouth with warm water. Then remove the excess gel with a toothbrush and toothpaste. The Home Bleaching products are simply cleaned with water.
7. You can apply the tooth whitening once a day until you have achieved the result you desire.

Recommendation

Before the first application, make sure you have no untreated caries or gingivitis. If unsure, consult a dentist before use. Please note that human teeth naturally differ in colour. A darker shade of tooth cannot be turned into a very light shade, not even after prolonged use. Take a comparative photo of your teeth before the first treatment to better assess the improvement.

If you experience a sensitivity reaction, such as irritated gums or pain-sensitive teeth, reduce the number and duration of treatments to reduce symptoms, without having to give up the treatment altogether. The sensitivity of gums and teeth is usually temporary and should subside completely within a few days after application.

Using special toothpaste for sensitive teeth also helps.

Please note that the use of the product is at your own risk and any liability for incorrect use will be rejected.

Warning

- Avoid contact with eyes
- Rinse eyes immediately if the product has come into contact with eyes.
- Consult your dentist before using for the first time
- Read internal user and safety information carefully before use.
- This product is not intended for use by persons under the age of 18.
- To be applied over a maximum 7 consecutive days and not more than once per day
- Occasionally irritated gums and pain-sensitive teeth may occur during application.
- Do not use in the case of untreated caries and known gingivitis.
- If you are allergic to one or more of the ingredients, please contact your dentist before use.
- Keep product away from children

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